

USER GUIDE

SUPER STROLLER® GUTTER WALKING FRAME

Code: 1846G



CONTENTS

1. Intended Use	2
2. Product Description	2
3. Definitions	2
4. Safety Warnings	3
5. Product Specifications	4
6. List of Components	4
7. Instructions for Use	5
8. Care and Maintenance	7
9. Limited Warranty	8

1. INTENDED USE

The folding gutter walking frame is designed for individuals who, due to a wide variety of reasons, require additional support while walking. It should be used as a tool to assist people who have mobility difficulties with standing, walking or sitting independently. It is designed exclusively for indoor use on solid, level surfaces and for persons with a maximum body weight of 100kg.

2. PRODUCT DESCRIPTION

This lightweight aluminium folding walking frame has large 6" wheels with a unique hand operated one touch brake system, padded seat, removable back rest, comes standard with basket, extra easy to fold for storage. Very easy to use in small areas.

3. DEFINITIONS

Symbols used in this user manual and their meanings.



WARNING! Failure to heed this warning may result in damage to the product or serious injury to the operator/user.



CAUTION! Beware of pinch points.



ATTENTION! Read and understand the instructions in the user guide before using the product.

4. SAFETY WARNINGS



ATTENTION! This user manual must be read before using the walking frame. Severe injury or death may result if user instructions, maintenance instructions and product warnings are not followed.



WARNING! The decision to use this walking frame should always be guided by a qualified Health Care Professional, who can also advise and assist with the correct height adjustment and usage instruction.



CAUTION! Fingers or other body parts can be trapped between the walking frame's folding mechanism during operation. Be extremely careful when adjusting the walking frame. Observe the dimensions of gaps between the walking frame's elements to avoid trapping your fingers or other body parts.

- DO NOT exceed the maximum weight capacity (including basket contents) or serious injury could result: Maximum weight capacity: 100kg
- DO NOT hang anything from the gutter frame other than the basket included. The basket has a weight limitation of 8kg
- Instruction and practice are necessary for successful use of the gutter frame.
- Ensure handbrakes are fully operational before each use.
- The gutter frame is to be used indoors on level, firm surfaces.
- Use the gutter frame only for the purpose it is intended for.
- DO NOT stand on the gutter frame.
- Ensure left and right handles are level to avoid instability problems.
- DO NOT use your gutter frame until the handles are properly adjusted and securely tightened.
- Always lock brakes when sitting on the gutter frame.
- DO NOT use the gutter frame as a wheelchair; doing so may cause a tip over, resulting in injury.
- DO NOT use the gutter frame to pull up on. This is not a transfer aid.
- DO NOT lean too heavily on the handles.
- DO NOT sit on the seat for long periods. (It is designed as a limited seating function for temporary rest)


5. PRODUCT SPECIFICATIONS

SPECIFIC DIMENSIONS

Maximum user weight	100kg
Approximate frame weight	7kg
Adjustable height – floor to trough	950mm - 1280mm
Width between handles	430mm
Seat height	550mm
Seat size	300mm x 300mm
Overall width	600mm
Inside width	440mm
Folded depth	320mm
Safe working load of basket	8kg
Usage	Indoor use only

PRODUCT LABELS

- Date of manufacture
- Model number

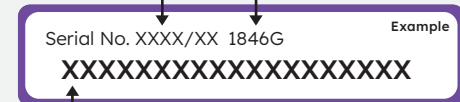

- Serial number

- Model name and number
- Refer to instructions for use



- Indoor use only
- Maximum safe working load
- Total width of walking frame



WARNING! All labels must remain readable. Immediately replace unreadable or missing labels. The serial number label is very important for identifying the product. It may not be removed.

6. LIST OF COMPONENTS

The walking frames are supplied in boxes and come folded, for easy transportation. The contents include one complete gutter frame folded, two gutters, one basket and one user manual.

ADDITIONAL OPTIONS



Pneumatic wheels (each)

Code: 5003FITTED



Single lever brake

Code: SLB



Oxygen bottle holder

Code: 1793



Folding tray (folds flat)

Code: 1847


7. INSTRUCTIONS FOR USE



CAUTION! Beware of pinch points

UNFOLDING YOUR WALKING FRAME FOR USE

Push down on the centre bar until frame opens (1) then put the basket in place (2). Fold down the seat (3). Ensure the walking frame is fully open (4).

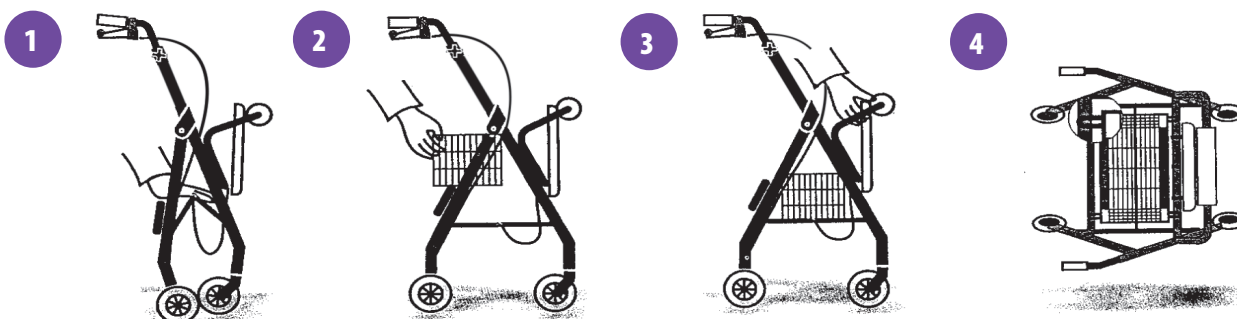


INSTALL BACKREST

Insert backrest into frame until the snap buttons engage. Ensure backrest is securely fastened

HOW TO FOLD

Lift the seat (5), remove the basket (6) and pull on the nylon rope (7) to fold the frame (8).

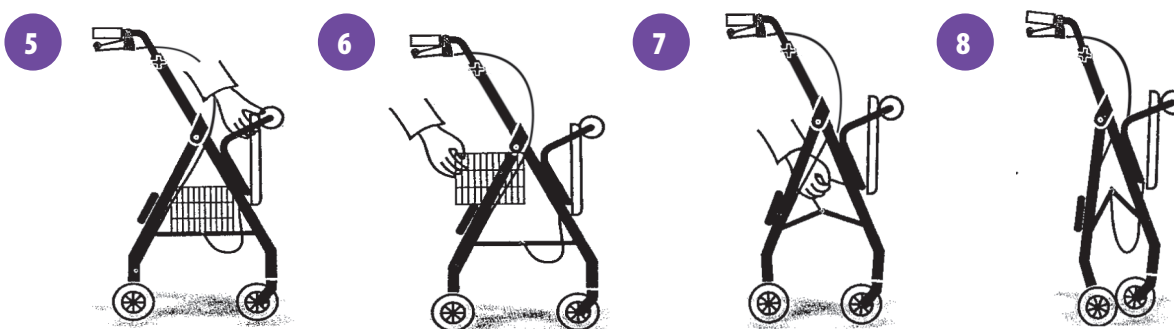



WARNING! Seek advice from a qualified health professional for the appropriate height/length level to suit your needs.



CAUTION! Beware of pinch points



STEM ADJUSTMENT

1. Remove adjustment knob and bolt from both stems. (9).
2. Insert stems into gutter frame (10).
3. Insert the bolt through the frame, stem at required height.
4. Ensure before tightening that the hexagonal bolt head is into the hexagonal hole in the frame.
5. Tighten knob so stems are firm.


FITTING GUTTER TO STEM

1. As per picture below fit the gutter component to the stem. (11).



USING THE HAND BRAKES

1. To operate the brake pull upwards on the handbrake lever (12).
2. When using the walking frame to move forward, the handbrake lever will remain in (13)
3. To lock the handbrake, push down on the handbrake lever until it "clicks" into the lock position (13) to (14).
4. To unlock the handbrake, pull up on the handbrake lever from (14) to (13).

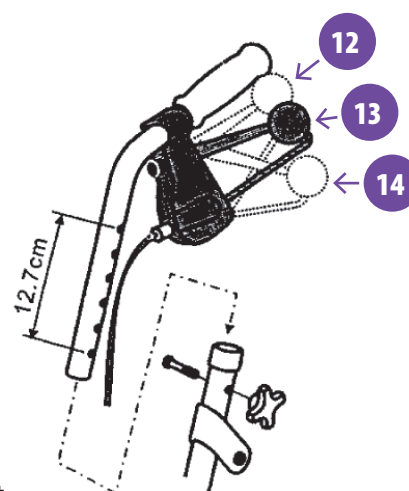

WARNING! When using the walking frame in a stationary position, the hand brakes **MUST** always be locked.

BRAKE TENSION ADJUSTMENT

1. To adjust cable at the handbrake handle, loosen lock nut and t cable adjuster in to loosen the brake or out to tighten the brake.
2. Bottom cable adjustment (above wheel)-turn nut towards brake to tighten or turn nut upwards to loosen brake.

8. CARE AND MAINTENANCE

CLEANING

- Clean the walking frame regularly with a damp cloth.
- Use a mild washing agent and warm water for excess dirt build up.
- Do not use harsh abrasives or bleach based agents which may cause damage or discolouration.

WHEELS

- Clean wheels with a mild washing agent and warm water.
- Lubricate with multi-purpose oil periodically or if wheels begin to squeak.

MAINTENANCE

- Ensure that every part of the walking frame is secure at all times.
- Check all screws for tightness. If not, tighten if necessary.
- Replace any broken, damaged or worn items immediately.

DISPOSAL

Products that can no longer be used are to be disposed of separately from household rubbish. This must be done according to the local and national regulations for environmental protection and raw material recycling. Please look at your local council's website for further information on the correct disposal.

MAINTENANCE SCHEDULE

Weekly	Periodically	Yearly	Serviceable Components
• Check brake function • Check all wheels, nuts and bolts are secure	• Check handgrips for wear • Check folding hinge bolts are secure • Check front and rear axle bolts are secure • Check frame for wear or cracks • Check wheels for damage or splitting • Check seat for wear and cracks. Ensure it is secured correctly to the frame	• Take your walking frame to an authorised dealer for a full check, clean and service.	• Hand grips • Wheels • Brake cables • Seat • Brake handles • Bearings • Nuts and bolts



WARNING! If faults found during any inspection remove from use immediately.



WARNING! Do not make any changes or modifications to this product without consultation from your Cubro agent.

9. LIMITED WARRANTY

Cubro Ltd warrants the product to be free from defects in materials and workmanship from the date of purchase for a period of 1 year (12 months). Please refer to cubro.co.nz for terms and conditions.

[illegible]

Super Stroller Gutter Walking Frame

USER GUIDE

[illegible]

[illegible]

© Copyright 2023 Cubro® Ltd

ALL RIGHTS RESERVED: No part of this publication shall be adapted, modified, reproduced, copied, or transmitted in any form or by any means including written electronic, mechanical, reprographic, photocopying, or recording means. Furthermore, this publication shall not be stored in whole, part, adapted, or modified form, in or for any retrieval system of any nature without the written permission of the copyright holder. Applications for authorisation of a reserved right of copyright owner shall be made in writing to the publisher. WARNING: The doing of any unauthorised act in relation to a copyright work may result in both a civil claim for damages and criminal prosecution.

09/06/23

Email sales@cubro.co.nz • Call us 0800 656 527 • Visit cubro.co.nz

Cubro®