

INTENDED USE

The ramp is intended to provide a gradual transition between surfaces of differing heights to assist individuals with reduced mobility. It is designed to facilitate access for mobility aids such as wheelchairs, walkers, and similar devices in indoor and outdoor environments.

SAFE USE



ATTENTION! This user manual must be read before using the ramp. Severe injury or death may result if user instructions, maintenance instructions and product warnings are not followed.

- **Never exceed the safe working load of 300kg.**
- Do not use if ramp is damaged or unstable.
- For wheelchairs and scooters only (not intended for commercial use).
- Ramp may be slippery when wet or icy.
- Surface temperatures can increase in direct sunlight.
- Before each use, check the ramp for worn, loose, or damaged parts. If you find any of these conditions, do not use the ramp. Contact your dealer for service or repairs, failure to do so may cause injury.
- Do not use the ramp if it is unstable.
- Do not tamper with, attempt to repair, or modify any portion of the ramp.
- Use caution when handling the ramp to avoid all potential pinch points.
- Do not use the ramp if its surface is covered with ice, snow, dirt, leaves, or other debris. Remove accumulations and sweep the surface clean before use.
- Always exercise caution when handling, installing, and or using the ramp.
- Use this ramp only for the purpose it is intended for. Improper use of the ramp could result in serious injury.



WARNING! Use only as instructed by your Health Care Professional.

PRODUCT SPECIFICATION

Code	557	557-900	557-1200	557-2000
Safe working load	300kg	300kg	300kg	300kg
Length	600mm	900mm	1200mm	2000mm
Side height	28mm	55mm	55mm	55mm
Inner side height	20mm	30mm	30mm	35mm
Fold side height	40mm	80mm	80mm	80mm
Fold width	400mm	380mm	380mm	370mm
Top lip height	20mm	30mm	30mm	35mm
Top lip length	85mm	100mm	100mm	100mm
Product weight	3kg	4.1kg	7.3kg	9.8kg

INSTRUCTIONS FOR USE

To open and then to use the ramp:


1. Stand the ramp on its end toward the ground.
2. Unfold the ramp in a similar manner to that of a book.
3. Whilst still supporting the ramp, stand to one side before gently lowering the ramp into position.

To stow / store the ramp:

1. Whilst standing to one side of the ramp, gently lift it into the upright position.
2. Close the ramp in a similar manner to that of a book.
3. Ideally store the ramp in its bag.



WARNING! Do not make any changes or modifications to this product without consultation from Cubro.

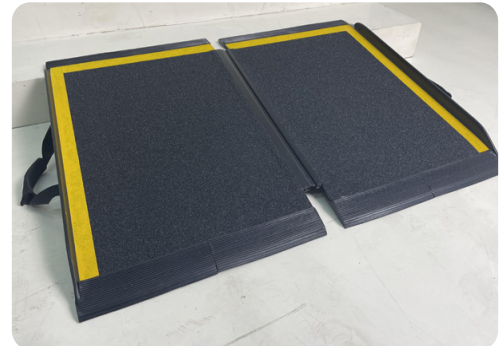

Cubro®
Strada



CARE AND MAINTENANCE

Inspection checks:

Before each use check the ramp for any obvious signs of damage. Particularly:

- a. The slip resistant surface.
- b. Dents / knocks / holes to the surface of the product that might cause a hazard.

Each month, stand the ramp on its end and:

- c. Check the hinge for any signs of significant wear or damage.
- d. Check that the slip resistant surface is still offering a degree of resistance to slip.
- e. Check that there are no knocks / dents / holes in the surface of the ramp that may cause a hazard.
- f. Check that all fixings are secure.
- g. Check that the handle.



WARNING! If faults found during any inspection remove from use immediately.

CLEANING

If the ramp needs to be cleaned, do so using only a mild detergent in warm water.

DO NOT use any cleaners/polishes that contain silicone.

DO NOT use any bleach or products containing bleach.

DISPOSAL

Products that can no longer be used are to be disposed of separately from household rubbish. This must be done according to the local and national regulations for environmental protection and raw material recycling. Please look at your local council's website for further information on the correct disposal.

SERVICING

For any servicing, maintenance and troubleshooting queries, please contact the Cubro Service Department:

New Zealand
0800 656 527 +64 7 578 7228
sales@cubro.co.nz
cubro.co.nz
149 Taurikura Dr, Tauriko, Tauranga NZ

Australia
1800 656 527
sales@cubro.com.au
cubro.com.au

For additional copies of this user guide or for a different format, please contact Cubro®.

LIMITED WARRANTY

Cubro Ltd warrants the product to be free from defects in materials and workmanship from the date of purchase for a period of 1 year (12 months). Please refer to cubro.co.nz for terms and conditions.

© Cubro Limited 2026

ALL RIGHTS RESERVED: No part of this publication shall be adapted, modified, reproduced, copied, or transmitted in any form or by any means including written electronic, mechanical, reprographic, photocopying, or recording means. Furthermore, this publication shall not be stored in whole, part, adapted, or modified form, in or for any retrieval system of any nature without the written permission of the copyright holder. Applications for authorisation of a reserved right of copyright owner shall be made in writing to the publisher. WARNING: The doing of any unauthorised act in relation to a copyright work may result in both a civil claim for damages and criminal prosecution.

21/04/2026